

Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

Specification Number:	Spec #004 Rev 06	Effective Date:	02-24-2025
FG Item Number:	30036	Allergens:	None
FG Item Description:	MIT45 Super K Extra Strong (30 mL) 12 Pack		

# MIT45 SUPER K EXTRA STRONG (30 ML) 12 PACK Finished Product Specification Sheet

Confidential

Batched Material Details		
Item Number:	10005	
Item Description:	Super K Extra Strong - Mixed Liquid	
Specification Number:	#B002	
Formula Number:	FOR10005	

Finished Product Details		
Product Type:	Dietary Supplement	
Storage Requirements:	Ambient Storage	
Shelf Life Info:	2 years past the first day of batching	
Additional Statement Claims:	N/A	

# Lot Code & Best By Date Code \*\*Image shown is for a different bottle; MIT45 (15 mL). The formatting is the Image: same, but the bottle color and size are different (the correct color is purple). The Lot Code is derived from the last 4 digits of the Batch Ticket Number For instance, if the Batch Ticket Number is MFG000123, then the product Lot Code shall be "0123". Sublots may be assigned in addition to the Lot Code per operational needs. Prefix sublots such as "00-0123" are determined by the distribution sales region. Suffix sublots such as "0123-1" should be Lot Code: assigned if there are significant changes with the components used or production process. Lot code may include both prefix and suffix "00-0123-1" if applicable. This shall be printed in yellow or white ink on the bottom of the bottle, and directly above the Best By Code. The Best By code for the product shall be 2 years past the first day of manufacture. For instance, if **Best By Code:** the batch was started on 01/01/2025, then the Best By Code shall be "01/2027". This shall be printed in yellow or white ink on the bottom of the bottle, and directly beneath the Lot Code.

MIT45.com 1 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

## Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

Bill of Materials					
Component	Item Number	Item Description Specification Qty		Qty per POS	UOM
Batch	10005	Super K Extra Strong - Mixed Liquid	#B002	360	mL
Bottle	01783	Painted (30 mL) Super K Extra Strong Bottle	1 #C080		ea
Lid	01784	Matte Black Dropper (30 mL) - CRC TE	#C081 12 e		ea
Unit Carton	01786	Diamond Box - Super K Extra Strong Individual unit (30 mL) #C083		12	ea
Label Sticker	03051	Super K Extra Strong Mitragynine Sticker	#C157 13		ea
Display Tray	01789	12ct Display Carton - Super K Extra Strong	#C086 1 e		ea
Tray Cover	01790	12ct Slip-Case Carton Cover - Super K Extra Strong (30 mL)	#C087 1 ea		ea
Wafer Seal	01382	Wafer Seal - Clear Acetate - 2 1/2" #C105 2		2	ea
Master Case	N/A	Master Case (19" x 13" x 10")	(19" x 13" x 10") N/A 0.083 e		ea
Таре	N/A	MIT45 Branded Tape	N/A	N/A	N/A

Compone	nt Images
Super K Extra Strong - Mixed Liquid	Matte Black Dropper (30 mL) - CRC TE
	의토임토외토웨

	Finished Product Manufacturing Instruction:
1	Supervise the bottles as they run through the line, ensuring that none fall over.
2	Check for loose or missing lids. Ensure that tamper resistant seal is intact.
3	Pack the bottles into their individual unit diamond boxes.
4	Place 12 of the diamond boxes into the display tray, all facing forward.
5	Place the tray cover over the filled tray, and seal with 2 clear round stickers on the bottom centers that are opposite of the notch on the assembled carton.
6	Pack the display boxes into master cases and tape with MIT45 branded tape. Place the appropriate sticker on the completed master case. Palletize according to specification.
7	When the pallet is complete, notify Warehouse Personnel and they will collect it.

MIT45.com 2 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

#### Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

Label Specifications		
Specification Number:	#C080	
Artwork Version:	V2	
QR Code Link:	MIT45.com (Home Page)	
Dimensions:	3.125" (H) x 1.32" (D)	
Material/ Finish:	Clear Glass Bottle / Matte Foil Finish	
- 4/-		

#### **Proof/Artwork Image**



#### WARNINGS

if adverse reaction occurs, stop use and consult a health professional. Do not use while pregnant or nursing. Consult with a health professional prior to use

Do not use with alcohol or prescription

Not intended for daily use.

These statements have not been evaluated by the US Food & Drug Administration. This product is not intended to diagnose, treat, cure, or prevent any disease or condition.



mit45.com 866-M/T-4555

380 W. Lawndale Dr. Salt Lake City, UT 84115



## Supplement Facts

Serving Size: 2 Droppers (2mL) Servings Per Container: About 15

#### **Amount Per Serving**

ratom Extract (Leaf)	188mg*
Mitragynine	25mg*
Speciogynine	<3.90%*
Speciociliatine	<6.50%*
Paynantheine	<7.90%*
7-OH Mitragynine	<1.00%*

Ingredients: Water, Glycerin, Kratom Extract, and Citric Acid Max Daily Servings: 8 (200mg).

V2

## **Supplement Facts Panel:**

#### Supplement Facts Serving Size: 2 Droppers (2mL) Servings Per Container: About 15 **Amount Per Serving** Kratom Extract (Leaf) 188mg\* Mitragynine 25mg\* <3.90%\* Speciogynine <6.50%\* Speciociliatine **Paynantheine** <7.90%\* 7-OH Mitragynine <1.00%\* Daily Value not established,

Ingredients: Water, Glycerin, Kratom Extract, and Citric Acid. Max Daily Servings: 8 (200mg).

#### WARNINGS

If adverse reaction occurs, stop use and consult a health professional.

Do not use while pregnant or nursing. Consult with a health professional prior to use.

Do not use with alcohol or prescription drugs.

Not intended for daily use.

These statements have not been evaluated by the US Food & Drug Administration. This product is not intended to diagnose, treat, cure, or prevent any disease or condition.

MIT45.com 3 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

## Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

	Unit Carton Specifica	ations	
Specification Number:	#C083		
Artwork Version:	V2		
UPC Code:	810149700071		
QR Code Link:	MIT45	5.com (Home Page)	
Dimensions:	4.25" (H	H) x 1.5" (W) x 1.5" (D)	
Material/ Finish:	Foil Board, Spot UV,	Emboss, Holographic Silver Paper	
	Proof/Artwork Ima	age	
53 <sub>G4</sub> Supplement Facts Supplement Fact	6.999	These statements have not been evaluated by the US Food & Drug Administration. This product is not intended to diagnose, treat, cure, or prevent any disease or condition.    Mit45.com 866-MIT-4555   380 W. Lawndale Dr. Salt Lake City, UT 84115   W. Lawndale Dr. Salt Lake City, UT 84	
upplement Facts Panel:	Supplement Facts Serving Size: 2 Droppers (2mL) Servings Per Container: About 15  Amount Per Serving  Kratom Extract (Leaf) 188mg*  Mitragynine 25mg* Speciogynine <3.90%* Speciociliatine <6.50%* Paynantheine <7.90%* 7-0H Mitragynine <1.00%*  *Daily Value not established.  Ingredients: Water, Glycerin, Kratom Extract, and Citric Acid. Max Daily Servings: 8 (200mg).	Warnings  If adverse reaction occurs, stop use and consult a health professional.  Do not use while pregnant or nursing.  Consult with a health professional prior to use.  Do not use with alcohol or prescription drugs.  Not intended for daily use.	

MIT45.com 4 of 13




Document Number: Form 007

Effective Date: 06-05-2024

Revision: 06

#### Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

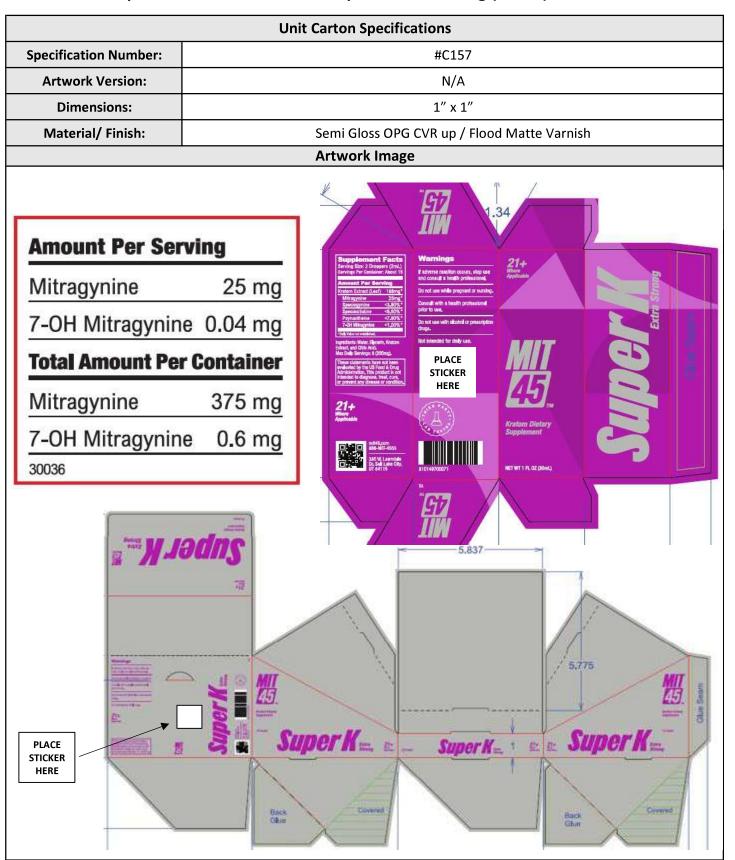




Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

#### Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

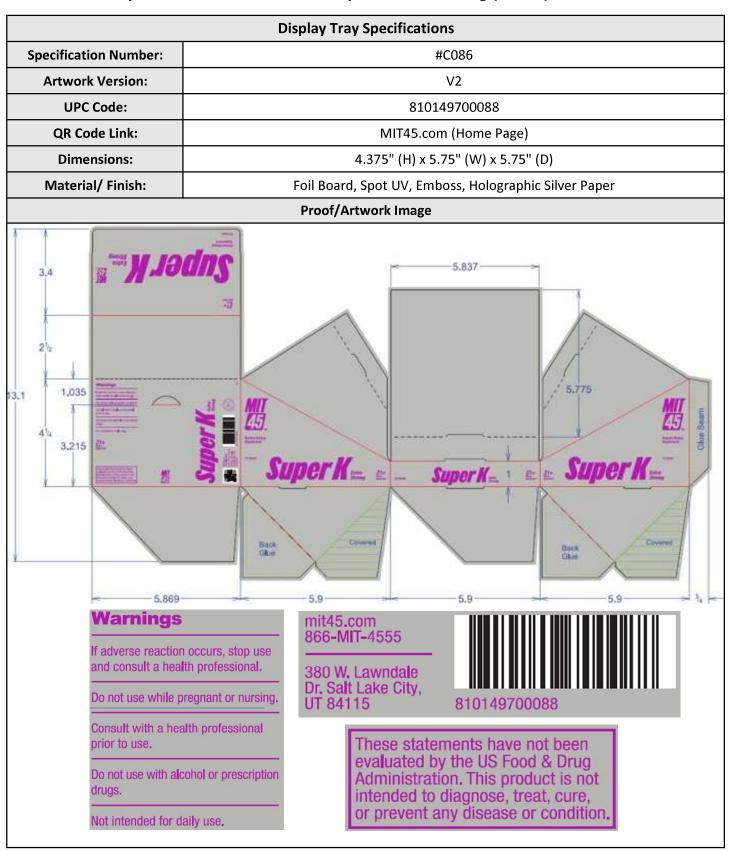

MIT45.com 6 of 13




Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

#### Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

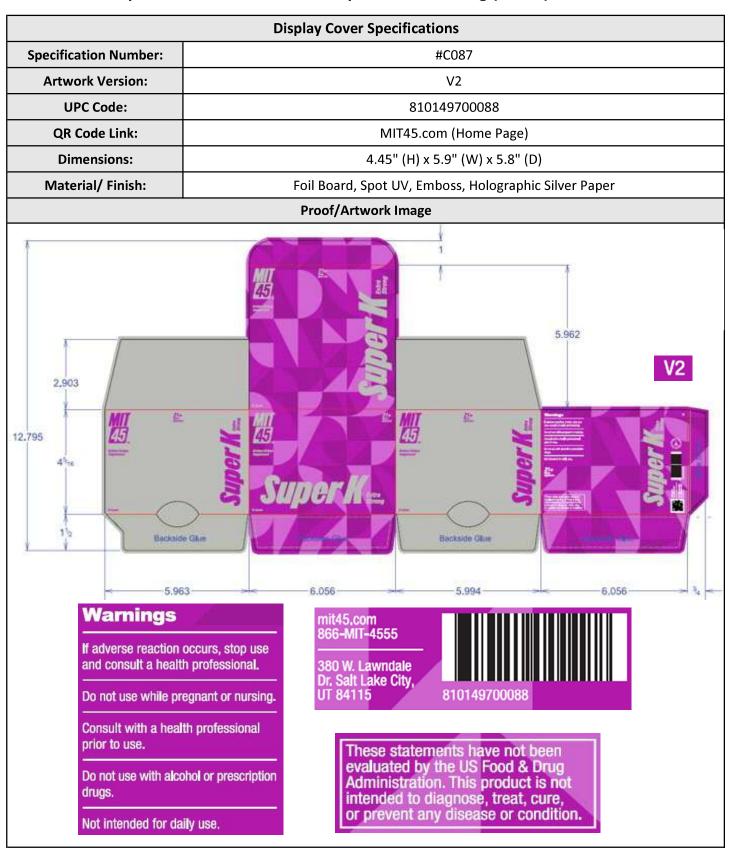

MIT45.com 7 of 13



**Document Number: Form 007** 

Effective Date: 10-01-2024

Revision: 07

## Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

Packout Specifications		
Packing Type	Quantity	Packout Description
		3 rows of 4 diamond boxes, all forward-facing. Seal the carton securely
Display Carton	12 bottles per carton	using two clear round stickers, placed at the center of the bottom on
		both the front and back sides (opposite of the notches).





Master Case (19" x 13" x 10")

12 cartons per master case

**2 layers of 6 display cartons,** with filler paper between the 2 layers, and taped with MIT45 branded tape.





**Pallet Configuration** 

30 boxes per pallet

**6 master cases per layer, 5 layers high.** Wrapped with minimum 3 layers stretch wrap and taped with MIT45 branded tape.



MIT MIT MIT MIT MIT MIT 45, 45, 45,

MIT45.com 8 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

## Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack



Weight Specs	
Minimum:	30 mL
Target:	31 mL
Maximum:	33 mL

MIT45.com 9 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

# Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

Acceptable Quality Limits
Critical (Zero Tolerance)
Empty containers
Leaking product
Incorrect component used
Critical (4%)
Illegible/Missing/Incorrect Lot Code or Shelf Life
Damaged components causing product exposure
Incorrect Label Artwork
Major (6.5%)
Damaged components creating esthetic defect
Print flaw/damage on the primary display panel (including primary artwork/SFP)
Packout count inaccuracy
Minor (10%)
Foreign matter or product residue on the outside of the component
Print flaw/damage on the non-primary display panel (excluding primary artwork/SFP)

Sample Requirements			
Type of Sample	Sample Quantity	Time on hand	
Testing Samples	3 bottles each pulled from the beginning, middle, and end of production (9 bottles total)	To be sent for testing as soon as possible after production completion.	
Product Retains	1 carton each pulled from the beginning, middle, and end of production (3 cartons total)	To be retained for no less than 1 year past the expiration / Best By date.	

MIT45.com 10 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

## Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

Testing Specifications				
Category	Name of Assay	Method	Limit	UOM
	ID Analysis	Cambium/FTIRR	≥ 95	%
Identity	Color/Appearance	Organoleptic (In-House)	Amber - Dark Green	N/A
	Scent	Organoleptic (In-House)	Herbal	N/A
	Composition	Organoleptic (In-House)	Liquid	N/A
Composition	Density (Specific Gravity)	Cambium/TM-017 ARL/2.006	1.13 – 1.17	g/mL
	Water Activity	Cambium/TM-017	Report Only	aw
	Total Aerobic Plate Count	Cambium/TM-060 ARL/USP<2021>	≤ 10,000	CFU/g
	Coliforms	Cambium/TM-059 ARL/AOAC 991.14	≤ 1,000	CFU/g
Purity -	Escherichia Coli (E. Coli)	Cambium/TM-059 ARL/USP<2022>	Not Detected	CFU/g
Pullty	Staphylococcus aureus	Cambium/TM-062 ARL/USP<2022>	Not Detected	CFU/g
	Salmonella	Cambium/TM-063 ARL/USP<2022>	Not Detected	CFU/g
	Yeast & Mold	Cambium/TM-061 ARL/USP<2021> or AOAC RI	≤ 1,000	CFU/g
	Mitragynine	Cambium/TM-052 UPLC-DAD	20 – 30	mg/serving
Strength -	Speciogynine	Cambium/TM-052 UPLC-DAD	≤ 3.90	%
Strength	Speciociliatine	Cambium/TM-052 UPLC-DAD	≤ 6.50	%
	Paynantheine	Cambium/TM-052 UPLC-DAD	≤ 7.90	%
Potential Contaminants	7-OH-Mitragynine	Cambium/TM-047 LC-TQ	≤ 1.00	%
Note	This is a natural product, variances may be found that are due to the growing and drying conditions, age, season, harvest time, geographic location, production process, etc.			

MIT45.com 11 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

## Spec #004 Rev 06 - MIT45 Super K Extra Strong (30 mL) 12 Pack

## **Review & Approval**

By signing below, I acknowledge that I have read, reviewed, understood, and approved of all the above specifications. I understand that only specifications from the most current signed revision shall be used for production of the above product, unless under a signed deviation.

Dept.	Title	Name	Signature	Date
lity	Document Control Specialist	Cherie Mayers	DocuSigned by: Cluric Mayers 87BC95A6A508472	2/20/2025
	Quality Supervisor	Lulu Brea	Signed by:  Ulu Bria  70AF16923CC1493	2/21/2025
Quality	Senior Quality Manager	Matthew Trimble	Signed by:  Matt Trimble  AD1BB8E654D946B	2/21/2025
	Vice President of Quality	Brian Anderson	Signed by: Brian Indurson 9FB5FE41C3DB488	2/20/2025
Operations	Production Supervisor	Fernando Vazquez or Blanca Perez	Signed by: Funando Vazquez C5A5A9124D49407	2/20/2025
Opera	Director of Manufacturing	Parker Brohl	Signed by:  Parker Brow  7497C538E91D403	2/20/2025
ERP	Senior Director of Operations	Jordan Briggs	Jordan Briggs C366AA73144A4B1	2/20/2025

MIT45.com 12 of 13



Document Number: Form 007

Effective Date: 10-01-2024

Revision: 07

Specification Number:	Spec #004 Rev 06	Effective Date:	02-24-2025
FG Item Number:	30036	Allergens:	None
FG Item Description:	MIT45 Super K Extra Strong (30 mL) 12 Pack		

# **Revision History**

Revision	Date	Description of Revision	Author
00	11-22-2022	New Specification	Courtney Burton
01	06-26-2023	Revised Form Template, Acumatica Stock Item # Update, Addition of Fill Weights, Update to Testing, & Addition of Pricing	Courtney Burton
02	12-12-2023	Update BOM descriptions and sample quantities	Courtney Burton
03	06-23-2024	Updated to new template. Correct item number for Display Carton. Added organoleptic image. Updated fill volume/weight limits. Removed corrective UPC sticker from packaging instructions. Updated Lot Coding to allow top of lid or bottom of bottle. Updated testing to correct density limits, remove Nickel, add Kratom sub-alkaloids, add Residual Solvents.	Cherie Mayers
04	11-18-2024	Updated entire specification for rebrand; all new artwork components, updated Lot Coding location, updated Strength testing limits (all alkaloids), added testing for pesticides, new packout instructions, and updated retention time.	Cherie Mayers
05	11-20-2024	Updated AQL for leaking product to be Zero Tolerance	Cherie Mayers
06	02-20-2025	Added Mitragynine stickers (CHG #25-006). Updated weight limits to volume limits (CHG #25-003), testing specifications (CHG #25-002), wafer seal sticker size (CHG #24-051), and packout specifications (CHG #25-005).	Cherie Mayers

MIT45.com 13 of 13